
Instructions for use  
Mode d'emploi  
Gebrauchsanweisung  
Instrucciones de uso  
Gebruiksaanwijzing  
Istruzioni per l'uso

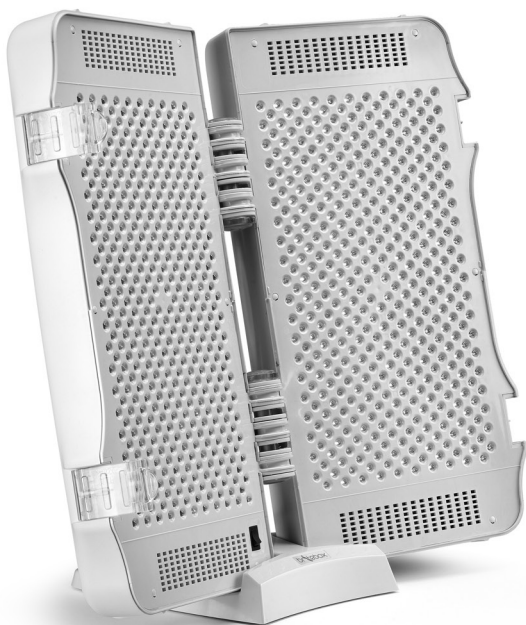English

Français

Deutsch

Español

Nederlandse

Italiano

# Languages

English.....	1
Français.....	10
Deutsch.....	19
Español.....	28
Nederlandse.....	37
Italliano.....	46

# Instructions for use

Please read these operating instructions before using the unit and save them for future reference.

## Introduction

Thank you for purchasing the britebox Light Therapy unit we wish you all the best with your treatments.

The britebox has been designed and manufactured specifically for the following treatments:

- Skin rejuvenation (anti-ageing) and Rosacea

Use britebox Red (633nm) and infrared (830nm)

Please note: Infrared is beyond human sight, however, by using any digital camera you will be able to see and check that the infrared LEDs are functioning correctly.

- Acne

Use britebox Red (633nm) and Blue (415nm)

Mary Read Ltd always appreciates customer feedback on our products and services so please do not hesitate to contact us... Thank you

## Recommended use

Make sure that the LED light source is approximately 8 inches (20cm) away from the area to be treated.

For all the above treatments it is recommended that you use the britebox for 15 minutes (any time of the day), for 5 days out of 7 and for two months. This is equal to 40 sessions of 15 minutes (total of 10 hrs).

The LEDs (panels) will start to lose their intensity after approximately 400 hours (Revive) and 200 hours (DermaClear). New panels can be bought from our website.

For skin rejuvenation and Rosacea we recommend that you have a top up treatment every 4 to 6 weeks and then repeat the 40 sessions each year.

For acne treatments you will need to monitor your condition and re-use when required.

There are no side effects to LED treatments.

The power consumption cost for each treatment is less than 1p.

## Results

### Results for Skin Rejuvenation and Rosacea

Results will be seen and felt within 4 – 6 weeks and will keep improving for 3 – 6 months after treatment has finished.

### Results for Acne treatment

A marked improvement will be seen within 4 -5 weeks with maximum results after 8 - 12 weeks. Please note that your condition may become slightly worse before improvements are seen.

# Desk Stand

For treatments using the britebox desk stand (enclosed with your britebox ). Simply place the desk stand onto a hard surface from which you wish the britebox to be used, with the 'britebox' label facing toward yourself (fig 1).

Unclip your britebox and open the two panels. Ensure that the switch is at the bottom of the unit (fig 2). Then place the britebox within the two slots within the desk stand and gently let the britebox slide into place (see fig 3).

fig 1

fig 2


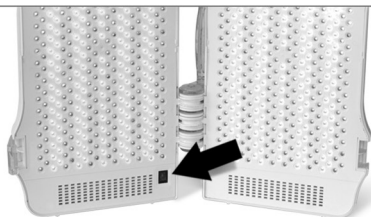
fig 3

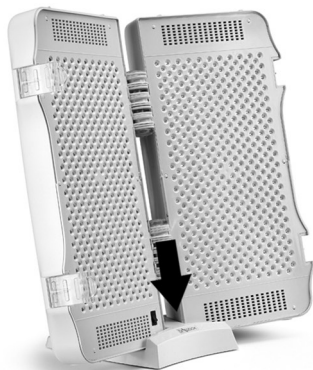Unpack your power supply and connecting lead (see fig 4). Plug the rounded barrel connector into the bottom side of the britebox (see fig 5). Plug in the IEC connector lead into the power supply and plug into the mains (see fig 6).

It is recommended that you exfoliate before each treatment for any of the above conditions.

Please note: It is essential that for acne treatments with the RED and BLUE LEDs that you use the goggles provided. This is due to the BLUE light (short wavelength) being detrimental to your eyes in extended use over time.

Goggles are not required for the RED and INFRARED treatments (unless preferred). However, you may find it more comfortable to keep your eyes closed during treatments.

fig 4

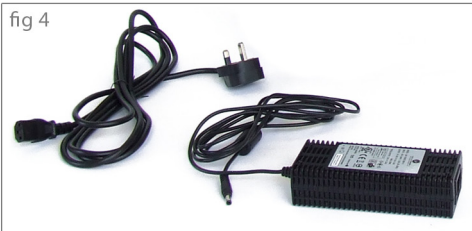
fig 5

fig 6

## OK, you are now ready for your first treatment...

Make sure that the LED light source is approximately 8 inches (20cm) away from the area to be treated.

Turn the britebox on by depressing the switch on the front of the unit. The LEDs will light, and your session begins. The britebox will automatically time 15 minutes and turn off. This marks the end of your treatment session. The fans will start when the britebox is up to temperature.

Relax for 15 minutes and let the britebox do the work...

After all treatments it is recommended that you moisturise.

## Maintenance

1. Always disconnect the unit from the mains when not in use.
2. Always unplug the britebox and let it cool down before any cleaning or maintenance. (No user serviceable parts inside).
3. Only clean the outside of the box with a damp cloth and use only a dry cloth to clean the LED panels.
4. For maintenance please contact Mary Read Ltd support department on +44 (0) 1949 83 83 83

Please note:

Do not use near water or within damp environments.  
Keep all vents and fan covers clear from any obstructions.  
Only use the goggles supplied.  
Only use the enclosed power supply

## Warranty & Guarantee

All Mary Read Ltd products come with a 60 day money back guarantee.

The britebox have a 2 year parts and labour warranty.

# Trouble Shooting

If your britebox is having technical difficulties then please follow the directions below.

Make sure that all leads are fully connected and the mains/britebox is switched on.

Check that the fuse has not blown in the mains plug.

If the britebox is receiving power (the LED power indicator will be on, (see fig 7.)) but the LED panels are not lit, then please refer to the error flash LED (top right panel (see fig 8)).

If the error flash LED is flashing then please note the number of flashes (see error flash table below).

If the error flash LED flashes 2, 3, or 4 times then please check if the vents are being obscured or, are you using the product within a high ambient/room temperature? (above 40degrees). If so then please use within an air conditioned environment or at a cooler time of the day

Error flashes on top right panel signal LED (see fig 8).

- 2 flashes Right panel over temperature (fan failure)
- 3 flashes Left panel over temperature (fan failure)
- 4 flashes Timer over temperature, fan failure or high ambient temperature
- 5 flashes Power supply over volts >25.5v
- 6 flashes Power supply under volts <21.5v or excessive load
- 7 flashes Internal error or excessive load

If the right panel signal LED flashes then please contact our customer support on: +44 (0) 1949 83 83 83.

fig 7


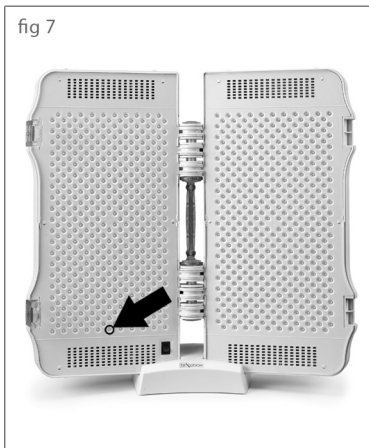
fig 8
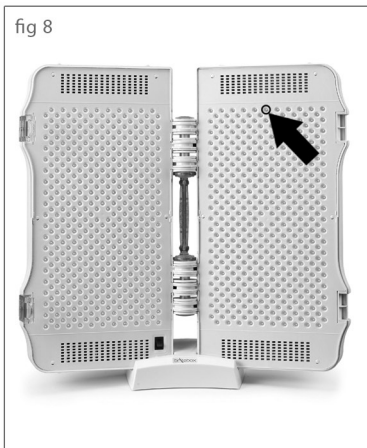
# Specifications

Mains Voltage 100 – 240V AC  
Max Current 2.5A  
Mains Plug Fuse Rating 5A  
Operating temperature -5°C to 35°C

No user serviceable parts

Please note: Only use supplied AC power supply.

## Light Output

Red – wavelength 633nm	(Power output 8.4mW/cm <sup>2</sup> )
Infrared – 830nm	(Power output 24mW/cm <sup>2</sup> )
Blue – 415nm	(Power output 14.6mW/cm <sup>2</sup> )

## Testing and Standards

Low Voltage Safety Assessment (LVD)

EN 60950-1: 2006

EMC

EN 61547:1995 AMD1: 2000

EN 55015:2006 AMD1: 2007

EN 61000-3-2: 2006

EN 61000-3-3: 1995 AMD1: 2001 AMD2: 2005

The britebox LED is not classed as a Medical Device.

# Customer Support

Mary Read Ltd, Tythby Grange, Bingham, Nottingham, NG13 8GQ

Phone Support: +44 (0) 1949 83 83 83

Email: [support@britebox-led.co.uk](mailto:support@britebox-led.co.uk)

## Copyright:

This material is copyrighted by Mary Read Ltd., and may be reproduced for internal use only. All other reproduction, in whole or in part, is prohibited without the written consent of Mary Read Ltd.

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© 2009 Mary Read Ltd. All Rights Reserved.
© 2009 Mary Read Ltd. Tutti i diritti riservati.